■ **GB** ELECTRIC TOOTHBRUSH

IE Operating instructions

FI SÄHKÖHAMMASHARJA

Käyttöohje

SE ELEKTRISK TANDBORSTE

FI Bruksanvisning

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NO ELEKTRISK TANNBØRSTE

Bruksanvisning

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
DE ELEKTRISCHE ZAHNBÜRSTE

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BALANCE

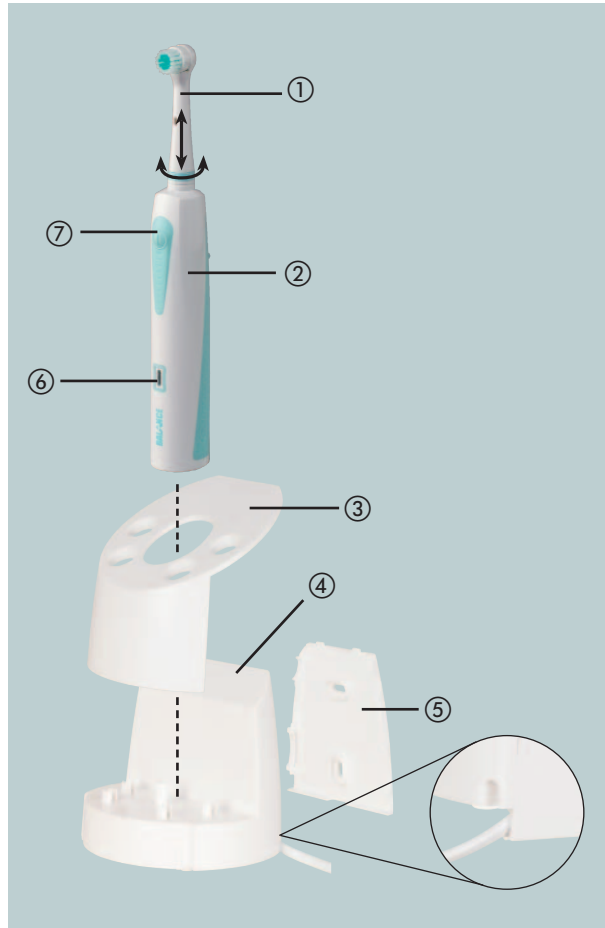Content	page
Safety instructions	4
Intended Use	4
Technical data	5
Included in delivery	5
Description of the appliance	5
Placement	6
Installing / removing brush heads	6
First use	7
Brushing one's teeth	7
Cleaning and Care	8
Disposal	9
Warranty and Service	9
Importeur	10

Read the operating instructions carefully before using the appliance for the first time and preserve this booklet for later reference. Hand over the manual to whomsoever you might pass the appliance on to.

Electric Toothbrush KH 107

Safety instructions

To avoid potentially fatal electric shocks:

- Ensure that the power cable can never become wet or damp during operation. Place the power cable such that it cannot become trapped or damaged in any way.
- Should the power cable become damaged, have it replaced by customer service before continuing to use the device.
- Never touch the charging station if it has fallen into water.
First of all remove the plug from the wall socket!
- Do not allow individuals (including children) to use the appliance whose physical, sensorial or mental abilities or lack of experience and knowledge hinder their safe usage of it if they have not previously been supervised or instructed in its correct usage.
Tell children about the appliance and its potential dangers so that they do not play with it.
- The electric toothbrush may only be recharged with the supplied charging station.
- You must under no circumstances immerse the hand appliance or the charging station in water or other liquids!

Intended Use

This appliance is intended for personal dental care in domestic situations. The appliance is not intended for use in commercial, industrial or medical environments. It is not suitable for animal care.

Technical data

Mains current:	220-240 V / 50 Hz
Power consumption :	1 W
Charging time:	app. 12 hours
Protection rating	
(Charging station):	IP X7
Protection class:	II 

Included in delivery

Electric Toothbrush

4 Brush heads

Storage box with fixture for two brush heads

2 Screws and 2 Plugs

Operating instructions


Description of the appliance

- ① Brush head
- ② Hand appliance
- ③ Brush holder
- ④ Charging station
- ⑤ Wall fixture
- ⑥ Display light
- ⑦ On/Off switch

Placement

The charging station ④ may only be operated free-standing in dry rooms. In moist rooms, especially in the vicinity of water (bath), the charging station ④ must be firmly fixed to the wall, to protect it from falling into water:

- Pull the wall fixture ⑤ downwards from underneath the charging station ④.
- Provide for the correct assembly of the wall fixture ⑤ close to a power socket. For this, use the plugs and screws supplied for this purpose.
- Subsequently, place the wall fixture ⑤, charging station ④ and brush holder ③ back together. Ensure that the power cable runs laterally through the recesses of the wall fixture ⑤ and the charging station ④. It may not become jammed!

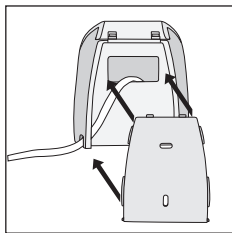Installing / removing brush heads

Turn the brush head ① clockwise and remove it upwards.

To install the brush head ①, apply the opposite procedure.

The brush head ① must click in on installation. Only use original brush heads, to avoid damage to the appliance.

(For replacement brush heads, see the order form)

First use

Before operating the appliance, ensure that it is in working condition. If not, it may not be used.

- Insert the power cord into the wall socket.
- Before the first use, clean the new brush heads ① with clear water.
- Place the hand appliance ② into the charging station ④ such that it easily slides into the guides in the plinth.
- Allow the hand appliance ② to charge up for ca. 24 hours before using it for the first time. After the first use, the battery requires only approximately 12 hours to completely recharge. The display light ⑥ on the front of the hand appliance ② glows red during the charging process.
- After every use place the hand appliance ② back into the charging station ④.

Brushing one's teeth

Warning: Due to the unaccustomed operation, some gum bleeding may initially occur, this is normal. Should this continue for more than two weeks, consult your dentist.

- Apply some toothpaste to the moistened brush.
- Switch the appliance on when it is in your mouth, to avoid splashing.
- To do this, press the On/Off switch ⑦.
- First of all brush the outer surfaces of your teeth. Place the brush head on the gums and direct the brush upwards or downwards over the tooth.
- Start again on the gums for each individual tooth.
- Clean the inner surfaces of your teeth with the same technique.

- To clean the chewing surfaces, direct the rotating brush head with light pressure over the chewing surfaces.
- Finally, thoroughly clean the back parts of the rear molar teeth.
- To switch off, press the On/Off switch again ⑦.

This appliance has a two minute timer. After two minutes the appliance becomes slower, indicating to you that the dental profession recommendation of two minutes has been complied with.

Cleaning and Care

- First remove the plug from the wall socket before cleaning the charging station ④.
- Never open the housing of the charging station or the tooth brush! These cannot be resealed in a water-tight manner.
- Only clean the charging station ④ with a soft, lightly moistened cloth. Never use detergents or solvents. These could damage the plastic surfaces.
- Clean the brush heads thoroughly under running water after use. We recommend cleaning the brush heads in a typical cleaning fluid for prosthetics every 4 - 8 weeks.
- The brush holder ③ can be cleaned without problem with water and some liquid soap. For this, the brush head holder can be simply removed by pulling upwards. However, allow all parts to dry completely before replacing them in the charging station ④.
- To maintain the batteries in the appliance, they should be completely charged up only when they become noticeably weaker - or at least every half-year.

Disposal

Do not dispose of the appliance in your normal domestic waste.

Dispose of the appliance through an approved disposal centre or at your community waste facility.

Observe the currently applicable regulations.

In case of doubt, please contact your waste disposal centre.

Dispose of all packaging materials in an environmentally friendly manner.

Warranty and Service

The warranty for this appliance is for 2 years from the date of purchase. The appliance has been manufactured with care and meticulously examined before delivery. Please retain your receipt as proof of purchase. In the case of a warranty claim, please make contact by telephone with our service department. Only in this way can a post-free despatch for your goods be assured.

The warranty covers only claims for material and manufacturing defects, not for wearing parts or for damage to fragile components, e.g. buttons or batteries. The appliance is intended solely for private domestic use, not for commercial applications.

The warranty is void in the case of abusive and improper handling, use of force and internal tampering not carried out by our authorized service branch. Your statutory rights are not restricted in any way by this warranty.

DES Ltd

Units 14-15

Bilston Industrial Estate

Oxford Street

Bilston

WV14 7EG

Tel.: 0870/787-6177

Fax: 0870/787-6168

e-mail: support.uk@kompernass.com

Kompernass Service Ireland

c/o Irish Connection

41 Harbour view

Howth Co. Dublin

Tel: 087-99 62 077

Fax: 18398056

e-mail: support.ie@kompernass.com

Importeur

KOMPERNASS GMBH

BURGSTRASSE 21

D-44867 BOCHUM

www.kompernass.com

Order card KH107

Quantity ordered	Name of article	Cost p. Unit	Total sum
	 Set of 8 Brushheads for: KH107 Toothbrush	£ 7,-	£

Payment method

Cheque payment

Creditcard/Debitcard

Bank Transfer

Add. £ 3,50 p&p.

Sender / Orderer

(please complete in block capitals)

Surname, first name

Street

Post code/town

Telephone

(Location, Date)

(Signature)

Order Card - How to order:

GB

1. Enter under "Sender/orderer" your name, address and telephone number (for possible queries) in block capitals.
2. Please tick your preferred payment method on the order card. You have got three options:

A: By cheque: Please place a cheque made payable to
DES Ltd

together with your completed order card into an envelope and post to our postal address below.

B: Creditcard/Debitcard for questions call:
0870 7876177

C: Bank Transfer (BACS): Please transfer the complete amount to our account below:

Barclays Bank

Windsor Branch

Account Number 20914959

Sort Code: 20-97-09

Domestic Electrical Solutions Ltd

On your transfer, please provide the article(s) you ordered as well as your complete name and address. Then post your completed order card to our postal address below.

Our postal address:

DES Ltd

Units 14-15

Bilston Industrial Estate

Oxford Street, Bilston

WV14 7EG

Important:

- Please affix sufficient postage to the letter.
- Write your name as the sender on the envelope.

Order card KH107

Quantity ordered	Name of article	Cost p. Unit	Total sum
	 Set of 8 Brushheads for: KH107 Toothbrush	8,00 €	€

€

Add. 2,00 € p&p.

Sender / Orderer

(please complete in block capitals)

Payment method

Cheque payment

Surname, first name

Street

Post code/town

(Location, Date)

(Signature)

Telephone

Order Card - How to order:

IE

1. Enter under "Sender/orderer" your name, address and telephone number (for possible queries) in block capitals. :

A: By cheque: Please place a cheque made payable to

Irish Connection

together with your completed order card into an envelope and post to our postal address below.

Our postal address:

Kompernass Service Ireland

c/o Irish Connection

41 Harbour view

Howth Co. Dublin

Important:

- Please affix sufficient postage to the letter.
- Write your name as the sender on the envelope.

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